

**File:** H460 293A exo mir21 1512.sds  
**Print Date:** Friday, December 15, 2023 15:35:26  
**User:** INSTR-USER  
**Plate Type:** Standard Curve  
**PCR Volume:** 10 µL

#### Document Information

**Operator:** INSTR-USER  
**Run Date:** Friday, December 15, 2023 12:06:58  
**Last Modified:** Friday, December 15, 2023 13:50:03  
**Instrument Type:** Applied Biosystems 7500 Fast Real-Time PCR System

#### Comments

SDS v1.5.1

#### Thermal Profile


Stage	Temperature	Time	Repeat	Ramp Rate	Auto Increment
1	50.0 °C	2:00	1	100	
2	95.0 °C	10:00	1	100	
3	95.0 °C	0:15	50	100	
	60.0 °C	1:00		100	

**Standard 7500 Mode**  
**Data Collection : Stage 3, Step 2**

#### Analysis Methods:

Delta Rn


#### Detector Information

Detector Name	Reporter	Quencher	Threshold	Baseline Start	Baseline End
 miR-21	SYBR	(none)	0.2	3	15

Well	Sample Name	Detector	Task	Ct	StdDev Ct	Quantity	Mean Qty	StdDev Qty	Filtered	Tm	User Defined #1	User Defined #2	User Defined #3
D2	H460 Exo 1	miR-21	Unknown	17.9005	0.124								
D3	H460 Exo 1	miR-21	Unknown	18.0238	0.124								
D4	H460 Exo 1	miR-21	Unknown	17.7764	0.124								
D5	H460 Exo 10	miR-21	Unknown	21.6017	0.154								
D6	H460 Exo 10	miR-21	Unknown	21.3444	0.154								
D7	H460 Exo 10	miR-21	Unknown	21.327	0.154								
D8	H460 Exo 100	miR-21	Unknown	Undet.									
D9	H460 Exo 100	miR-21	Unknown	Undet.									
D10	H460 Exo 100	miR-21	Unknown	Undet.									
E2	H460 Exo 1K	miR-21	Unknown	Undet.									
E3	H460 Exo 1K	miR-21	Unknown	Undet.									
E4	H460 Exo 1K	miR-21	Unknown	Undet.									
E5	H460 Exo 10K	miR-21	Unknown	Undet.									
E6	H460 Exo 10K	miR-21	Unknown	Undet.									
E7	H460 Exo 10K	miR-21	Unknown	Undet.									
F2	293A Exo 1	miR-21	Unknown	21.417	0.0234								
F3	293A Exo 1	miR-21	Unknown	21.4093	0.0234								
F4	293A Exo 1	miR-21	Unknown	21.3731	0.0234								
F5	293A Exo 10	miR-21	Unknown	25.2716	0.195								
F6	293A Exo 10	miR-21	Unknown	25.1876	0.195								
F7	293A Exo 10	miR-21	Unknown	24.8994	0.195								
F8	293A Exo 100	miR-21	Unknown	Undet.									
F9	293A Exo 100	miR-21	Unknown	Undet.									
F10	293A Exo 100	miR-21	Unknown	Undet.									
G2	293A Exo 1K	miR-21	Unknown	Undet.									
G3	293A Exo 1K	miR-21	Unknown	Undet.									
G4	293A Exo 1K	miR-21	Unknown	Undet.									

G5	293A Exo 10K	miR-21	Unknown	Undet.
G6	293A Exo 10K	miR-21	Unknown	Undet.
G7	293A Exo 10K	miR-21	Unknown	Undet.


## Raw Data

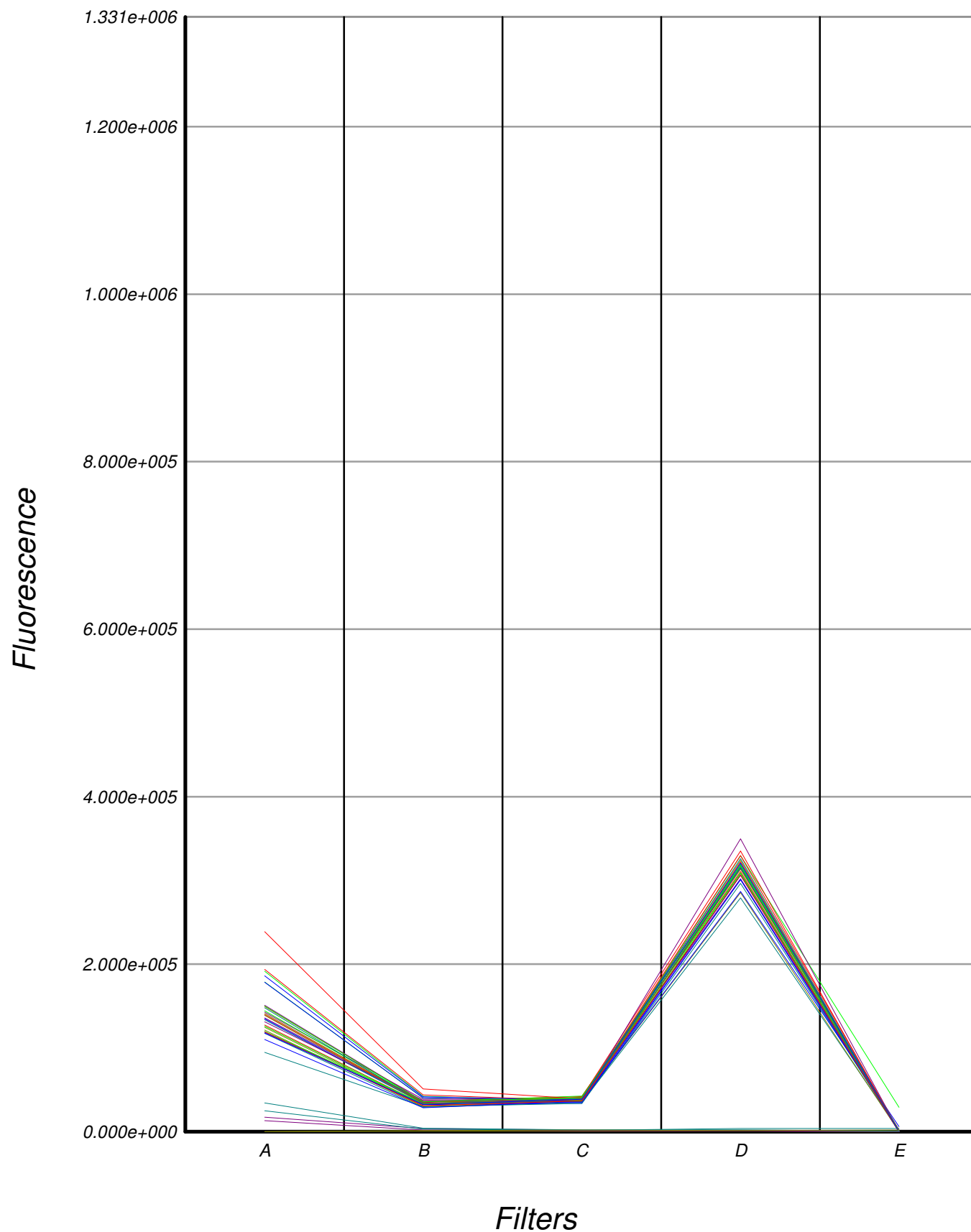

Reading: 1

Well(s): D2-D10,E2-E10,F2-F10,G2-G10

## Delta Rn vs Cycle

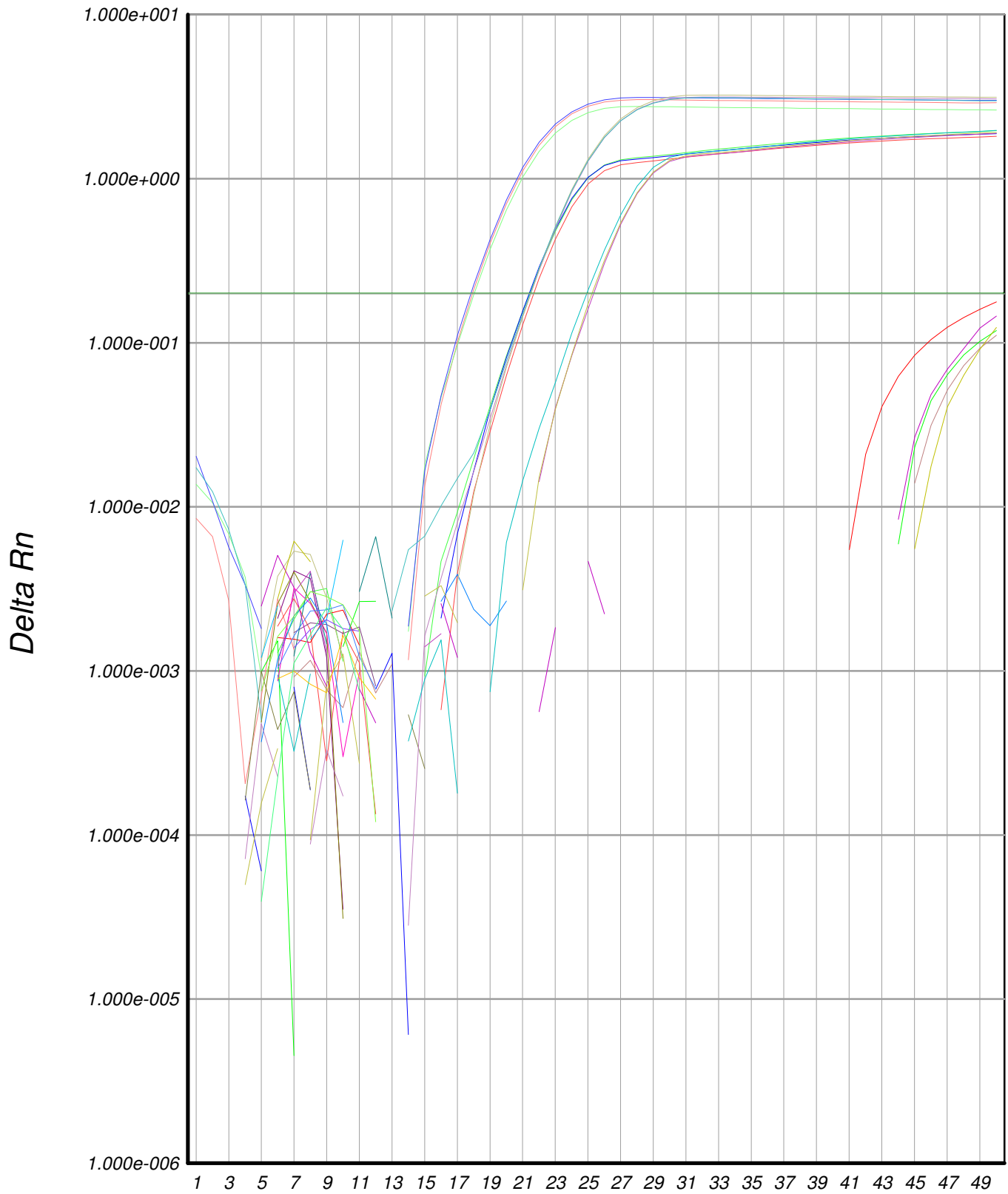

Selected Detector: All

Well(s): D2-D10,E2-E10,F2-F10,G2-G10